

PATIENT EDUCATION

Usage of Syringe for G-Tubing Venting

7-01-043



Properly administering your medication at home is important to your safety. In the event of an emergency, call 911.

SUPPLIES:

- 60ml syringe
- "Y" connector

PROCEDURE:

1. Clean work area. Wash hands thoroughly for at least 20 seconds. Gather supplies.
2. Place a "Y" connector into the feeding tube.
3. Place the tubing from the feeding bag into one side of the "Y."
4. Place a 60ml catheter tip syringe (without the plunger) into the other side of the Y. Venting syringe must be kept elevated 4 – 6 inches above the stomach to avoid a siphon effect. Wrap tape around the top of the syringe allowing enough tape to make a tab that can be pinned to the infants clothing or suspended from the crib rail.
5. Some patients will be ordered to remain vented at all times, while others may just be vented after a feeding then clamped. Follow the orders provided from you doctor.

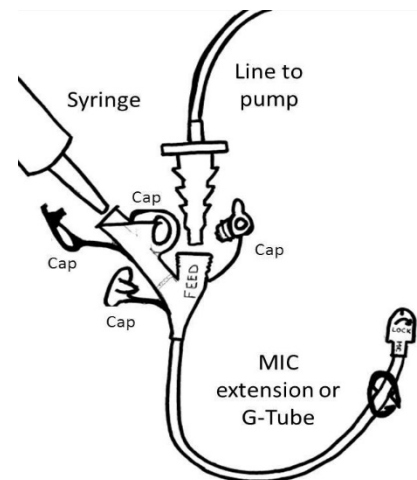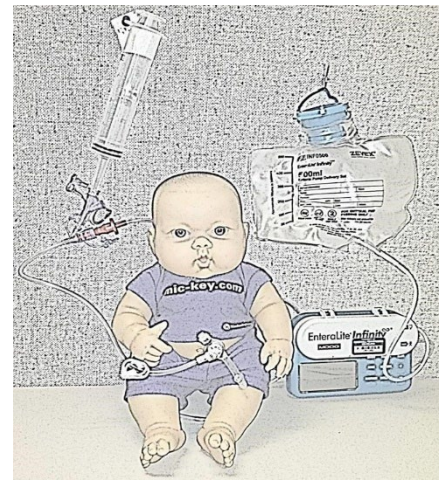

CAREPATHrx

The information included in this document is intended for informational and educational purposes only. It is not intended to replace medical advice offered by a physician or other qualified health provider. If you think you are experiencing a medical emergency, call 911 immediately.

Copyright© 2021 by CarepathRx. All rights reserved.